



Non-Adherent Dressings




# Pain-free removal.

## Telfa Non-Adherent Dressings

- Feature a perforated polyester film, containing a dense, non-woven cotton absorbent core that can be used on a variety of wounds
- Have a non-adherent surface with smooth aperture edges for intact and pain-free removal
- Allow exudate to pass through optimal size apertures without allowing tissue to adhere to the dressing
- Will not adhere to sutures, staples or adhesive wound closures


### Under the Microscope

Here you can see how Telfa dressing's superior apertures compare to our competitors.<sup>1</sup> Notice the optimal size apertures of the Telfa dressing, which allow exudate to pass through without allowing tissue to adhere to the dressing.

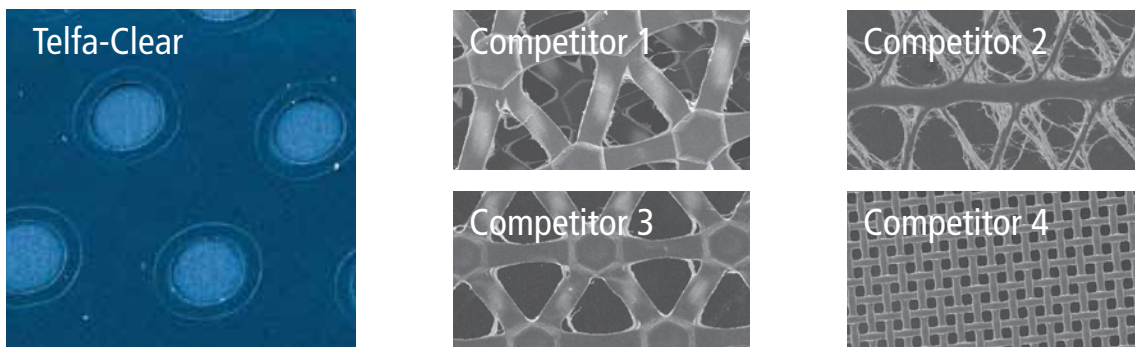

### A Comparison of Adherency Among Dry, Non-Adherent Dressings


In this test, adherence was measured by the amount of force required to remove the dressing from a two-inch strip of Bovine liver.<sup>1</sup>

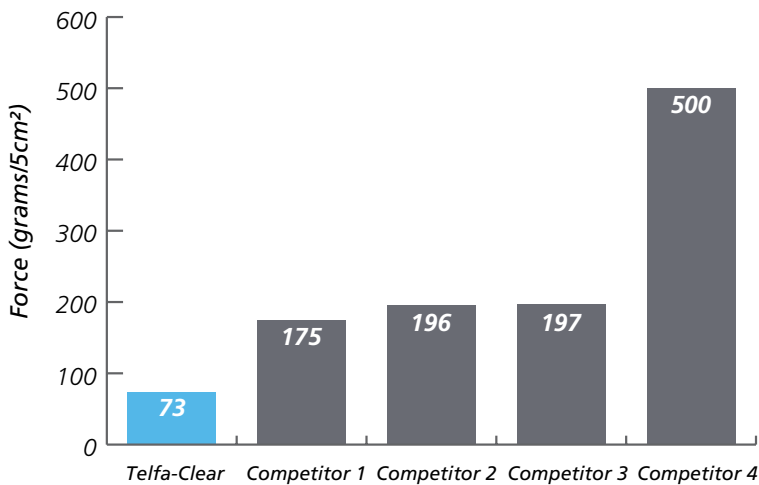

## Telfa Clear Dressings

Non-Adherent perforated polyester film


|                      |  |
|----------------------|--|
| <b>Non-Adherent</b>  | Allows for new tissue growth without interrupting the healing process. Allows for pain-free removal; Will not adhere to sutures, staples, or adhesive wound closures |
| <b>Non-Absorbent</b> | When combined with an appropriate secondary dressing, may be used on lightly, moderately, or heavily draining wounds   |
| <b>Clear Film</b>    | Allows assessment of the wound site without disrupting the healing process; May be left in place for several days  |
| <b>Versatile</b>     | Available in many convenient sizes; May be used on a variety of wounds   |
| <b>May be Cut</b>    | Allows for custom wound applications   |

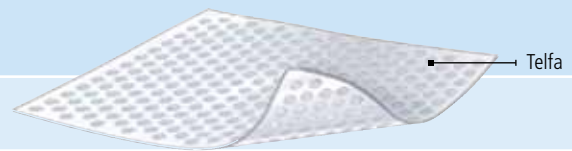

## Telfa "Ouchless" Non-Adherent Dressings

Non-adherent perforated polyester film, plus absorbent pad


|                       |  |
|-----------------------|--|
| <b>Non-Adherent</b>   | Allows for new tissue growth without interrupting the healing process; Allows for pain-free removal; Will not adhere to sutures, staples, or adhesive wound closures |
| <b>Absorbent Core</b> | Cotton core offers light to moderate absorbency to remove excess fluid from wound site   |
| <b>Versatile</b>      | Available in many convenient sizes; May be used on a variety of wounds   |

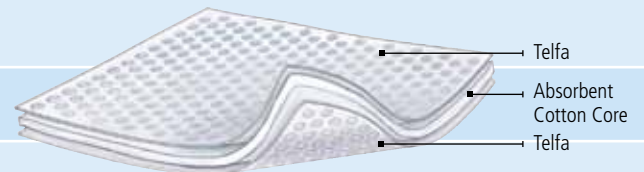

## Telfa Island Dressings and Telfa Plus Island Dressings

Non-Adherent perforated polyester film, plus pad with adhesive border

|                               |  |
|-------------------------------|--|
| <b>Non-Adherent</b>           | Allows for new tissue growth without interrupting the healing process; Allows for pain-free removal; Will not adhere to sutures, staples, or adhesive wound closures |
| <b>Absorbent Core</b>         | Cotton core offers light to moderate absorbency to remove excess fluid from wound site (Telfa Plus offers thicker pad for greater absorbency)                        |
| <b>Strike-Through Barrier</b> | Unique adhesive border helps prevent fluid strike-through and resists external bacteria penetration  |
| <b>Soft Non-woven Backing</b> | Conformable; comfortable to wear   |
| <b>Adhesive Border</b>        | No requirement for secondary dressing for fixation; Saves nurses' time. Ideal dressing for surgical wounds   |
| <b>Versatile</b>              | Available in many convenient sizes; May be used on a variety of wounds   |

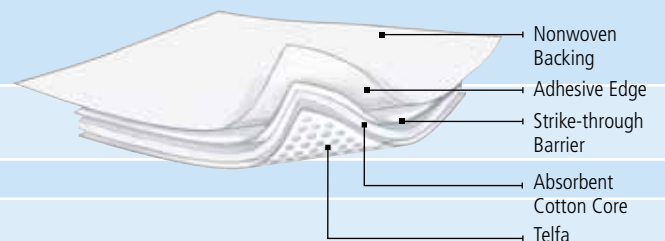

### Telfa Clear Wound Dressings

| Code                             | Description | Per Box | Ship Case |
|----------------------------------|-------------|---------|-----------|
| Sterile 1's in Peel-Back Package |             |         |           |
| 1109                             | 3" x 3"     | 25      | 50        |
| 1111                             | 4" x 5"     | 25      | 50        |
| 1113                             | 12" x 12"   | 25      | 50        |
| 1114                             | 12" x 24"   | 25      | 50        |

### Telfa Clear Wound Dressings (Bulk)

| Code            | Description   | Ship Case |
|-----------------|---------------|-----------|
| Nonsterile Bulk |               |           |
| 1115            | 39" x 25 yds. | 4 rolls   |

### Telfa "Ouchless" Non-Adherent Dressings

| Code                             | Description | Per Box | Ship Case |
|----------------------------------|-------------|---------|-----------|
| Sterile 1's in Peel-Back Package |             |         |           |
| 1050                             | 3" x 4"     | 50      | 900       |
| 1169                             | 3" x 6"     | 50      | 750       |
| 1238                             | 3" x 8"     | 50      | 600       |
| 28959                            | 8" x 10"    | 25      | 400       |

### Sterile 1's in Peel-Back Package, Perforated Carton for Easy Dispensing

|      |         |     |       |
|------|---------|-----|-------|
| 1961 | 2" x 3" | 100 | 2,400 |
| 2132 | 3" x 4" | 100 | 2,400 |

### Nonsterile Bulk, Non-Adherent Film on One Side

|      |          |     |       |
|------|----------|-----|-------|
| 2891 | 3" x 8"  | 200 | 1,000 |
| 3279 | 8" x 10" | 125 | 500   |

### Nonsterile Bulk, Non-Adherent Film on Both Sides

|       |         |  |       |
|-------|---------|--|-------|
| 12388 | 3" x 8" |  | 4,000 |
|-------|---------|--|-------|

### Telfa Sterile Strips

| Code                   | Description | Per Box | Ship Case |
|------------------------|-------------|---------|-----------|
| Sterile Telfa Dressing |             |         |           |
| 28959                  | 8" x 10"    | 25      | 400       |

### Telfa Island Dressings

| Code                             | Description | Per Box | Ship Case |
|----------------------------------|-------------|---------|-----------|
| Sterile 1's in Peel-Back Package |             |         |           |
| 7539LF                           | 2" x 3¾"    | 50      | 400       |
| 7540                             | 4" x 5"     | 25      | 200       |
| 7541                             | 4" x 8"     | 25      | 100       |
| 7542                             | 4" x 10"    | 25      | 100       |
| 7544                             | 4" x 14"    | 25      | 50        |
| 7550                             | 4" x 4"     | 25      | 200       |
| 7551                             | 6" x 6"     | 25      | 100       |

### Telfa Plus Island Dressings

| Code                             | Description | Per Box | Ship Case |
|----------------------------------|-------------|---------|-----------|
| Sterile 1's in Peel-Back Package |             |         |           |
| 2562                             | 4" x 6"     | 25      | 100       |
| 2563                             | 6" x 7"     | 25      | 100       |
| 2564                             | 6" x 10"    | 25      | 100       |
| 2565                             | 8" x 8"     | 15      | 60        |

### Telfa "Ouchless" Adhesive Dressings

| Code                             | Description | Per Box | Ship Case |
|----------------------------------|-------------|---------|-----------|
| Sterile 1's in Peel-Back Package |             |         |           |
| 6017                             | 2" x 3"     | 100     | 2,400     |
| 7643                             | 3" x 4"     | 100     | 1,200     |

<sup>1</sup> Quigg, J et al, Technical and Clinical Experience with Telfa Clear, 1994, Kendall Healthcare Products Co.

COVIDIEN, COVIDIEN with logo, Covidien logo and positive results for life are U.S. and internationally registered trademarks of Covidien AG. Other brands are trademarks of a Covidien company. © 2011 Covidien.

Manufactured by: Covidien Inc, 15 Hampshire Street, Mansfield, MA, USA, 02048

CDN-P00175-E (H6313) Rev. 2011/11

